smilesonic **KIDS**

USER MANUAL

CHAPTER 1

Meet your new friends

Hi, it's nice to meet you!

Starting today, we will take care of your teeth together.

First, pay attention to the elements of your toothbrush.

Delicate bristles that ensure sensitive and pleasant cleaning

Replaceable head

Handle

Power button

Battery and mode indicator
(In the Clean mode, the light is white.
In the Sensitive mode, the light is blue.
The light flashes when there is a need to charge the device.)

Stand with the induction charger

CHAPTER 2


Teeth love to be brushed

Let's prepare for brushing

1

Soak the head
in water

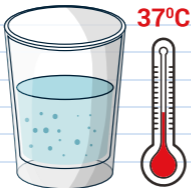2

Prepare a glass of
warm water


3

Put a bit of toothpaste
on the head

It's time to start brushing!

Put the head of the toothbrush in your mouth and press the power button to switch it on. Before you start brushing, select the mode appropriate for you using the power button (Clean – standard movement speed; Sensitive – more delicate movement speed).

The toothbrush will measure the time of brushing and will switch off automatically.

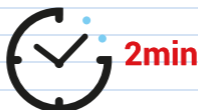Is it that simple?

Separate your teeth into 4 equal sections, and brush each section for 30 seconds. Slowly and delicately, move from one tooth to another with your toothbrush. After brushing, thoroughly rinse your mouth and teeth. Don't forget to clean your toothbrush and rinse the head.

We work towards healthy gums

If you are just beginning to use a sonic toothbrush, brush your teeth using the Sensitive mode. Don't press your toothbrush against your teeth too hard and avoid contact with gums. After brushing, place the toothbrush somewhere safe.

PSST! LET ME TELL YOU A SECRET!
YOUR TEETH WILL BE HAPPIEST
WHEN YOU BRUSH THEM
EVERY MORNING AND
EVERY EVENING.

CHAPTER 3

The things our teeth enjoy apart from brushing.

1. Brush your teeth every morning and every evening

Regular brushing yields the best results. Invest in white teeth and fresh breath.

2. Choose a head that was designed exclusively for children

The teeth and gums of children are delicate and sensitive. Select a head with appropriate bristles.

3. Select the appropriate mode

If you are just beginning to use an electric toothbrush, if your teeth are particularly sensitive (pain or bleeding appears when brushing), or if you simply prefer more delicate cleaning, use the Sensitive mode. In all other cases, use the Clean mode.

4. Use toothpaste for children

The teeth of children differ from those of adults. Use special toothpaste for children.

5. Teeth prefer delicate brushing

Don't press the toothbrush against your teeth too hard and avoid contact with gums.

6. Rinse your mouth with water after meals

After eating, the enamel on your teeth may be weakened. After you finish a meal, rinse your mouth and wait 20 minutes before brushing.

7. Drink more water

More humidity in the mouth means better protection against bacteria.

CHAPTER 4

Answer me! Children ask questions

Question: When can I start brushing my teeth with an electric toothbrush?

Answer: When we are 2-3 years old we already have about 20 milk teeth. A child that is 3 years old can already be given a gift in the form of an electric toothbrush. At that time, the child should also be taught further healthy habits improving the hygiene of the oral cavity.

Question: Why is my toothbrush for children different than my mom and dad's?

Answer: Teeth and gums are more sensitive in children than in adults. The speed of the toothbrush for children is lower, and the bristles of the head are soft and sensitive.

Question: Why do I have to brush my teeth? My milk teeth will be gone anyway.

Answer: Everything starts with milk teeth. It is at the moment when they start to appear that the good habits – prophylaxis, hygiene and a diet – start to form in children and adults. They are of great importance for the health of children in later years.

Question: Why we shouldn't brush our teeth immediately after a meal?

Answer: After a meal, we should wait for about 20 minutes before brushing our teeth. It is important to remember that food often includes acids, which soften and weaken the enamel. After eating, it is enough to rinse our mouth with cold water.

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Technical details

Product name: **Smilesonic Kids**

Brushing modes: **2 modes – Clean (62000 brush movements per minute) and Sensitive (52000 brush movements per minute)**

Operating voltage: **4 V**

Power rating: **3 W**

Total working time: **60-180 minutes**

Charging time: **6-15 hours**

Battery: **lithium battery/14500/650 mAh**

Waterproof rating: **IPX7**

The set contains: toothbrush handle, 2 heads, charger

FAQ

1. Why does my toothbrush vibrate slowly or the vibrations get weaker?

Answer: The head of the toothbrush might be incorrectly or too tightly attached to the handle. Detach the head and attach it again.

2. Why does my toothbrush lose its charge so quickly after full recharging?

Answer: Make sure the toothbrush is positioned vertically when charging. If the problem persists, contact us to repair or replace the device.

3. Why doesn't my toothbrush switch on?

Answer: Make sure that the toothbrush was properly charged and continue using the device normally after charging. If problems persist, contact us to repair or replace the device.

4. Is it possible to replace the power cord?

Answer: No, if the cord is damaged, you have to stop using the device and dispose of it.

5. Is the charger waterproof?

Answer: The charger features IPX7 waterproof class, but it is important not to keep it soaked or submerged in water. Make sure the charger is dry before connecting it to the power outlet and initialising recharge.

Proper brushing techniques

1. Slowly and gently move the brush head along the teeth in a back-and-forth motion, allowing the bristles to get to the interdental surfaces.
2. Do not brush the teeth the same way you did while using the manual toothbrush.
3. In order to increase the effectiveness of brushing, gently apply the toothbrush to the teeth, allowing it to brush them for you.
4. It is recommended to brush not only your teeth, but gums as well. The suggested order of brushing is: the outer surface, the inner surface and finally the chewing surface.

Cleaning the toothbrush

1. Clean the toothbrush regularly. Wash the handle and the head with running water after each use.
2. Carefully dry all parts before putting them on the stand.
3. Before cleaning the charger, make sure to unplug it
4. It is forbidden to clean the charger under running water or soaking it.
5. Carefully dry the charger before plugging it back on

Charging

1. The toothbrush included in the set is charged and ready for use. Please use it until a full discharge.
2. Rapid blinking of the battery indicator means that the battery power is running low.
3. Plug in the charger to mains power supply and then place the toothbrush on it.
4. During the charging procedure, the battery indicator should blink.
5. The indicator will continue blinking until the toothbrush is fully charged.
6. A full charging cycle lasts from 6 to 15 hours, allowing for 60 to 180 minutes of continuous work (depending on the selected brushing mode).
7. Once charged, the device will indicate the full battery by a constant light.
8. We can unplug the charger or leave it be (the device is protected against overcharging).

Safety

1. Before using the device for the very first time, make sure to read its instruction manual.
2. Before plugging the device, make sure that the voltage shown on the device is the same as the one in the outlet.
3. If the device is damaged in any way or makes unusual noises after being powered on (head, handle, charger, travel case), do not use it. It is forbidden to try to repair it on your own. The device does not contain any elements made for the user's repairs.
4. If the charger, travel case or brush heads are damaged, replace them with the original parts of the same type in order to avoid any potential danger.
5. While unplugging the device, always grab and pull the plug, not the cable.
6. The charger and power cord should be used and stored away from heat sources, water and damp.
7. The device should not be used by people (including children) with physical, sensory or mental impairments and those lacking knowledge or experience necessary to operate such devices. Children should be instructed not to use the device as a toy.
8. The device should not be used in any way other than the one specified in the manual. The tooth brush is made only for cleaning teeth, gums and tongue.
9. The manufacturer does not hold responsibility for damage caused by improper use of the product.
10. Do not use unofficial accessories.
11. Do not use the toothbrush and the charger while in bath or when taking a shower.
12. Do not use this device in places where it can be placed in or thrown into water.
13. Do not submerge the handle or the charger in water or any other liquid. After cleaning the device, make sure that every component is dry before plugging it back.
14. The device should not be used with wet hands.

15. The device is for home use only. It is not made for professional use.
16. If you had to undergo a jaw/gum surgery within the last 2 months, make sure to contact your doctor before using the device.
17. Consult your dentist if your gums start bleeding after using the toothbrush.
18. If you are using a pacemaker or a similar device, contact your doctor or the product's manufacture before using the toothbrush.
19. Replace the brush head every 3 months or sooner if you notice signs of wearing out.
20. Do not clean the brush heads in the dishwasher.
21. Do not use heads other than the ones recommended by the manufacturer.

**Correct disposal of this product
(waste Electrical & electronic equipment)**

After the product has reached the end of its service life, it is forbidden to throw it away along with regular waste. Instead, it should be brought to your local collection/recycling point for electric and electronic devices. This procedure is indicated by the symbol found on the product, packaging or in the manual. The materials used for manufacturing the device can be recycled and used again. Thanks to recycling, or other forms of disposing the device, you are all contributing to protecting the environment. The information about finding the correct recycling point can be obtained at your local municipal administration or the retail shop where you purchase the device.

smilesonic Ltd

C220d Trident Business Centre
89 Bickersteth Road
SW17 9SH London
United Kingdom

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