# USER MANUAL

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### **smilesonic Ltd**

C220d Trident Business Centre  
89 Bickersteth Road  
SW17 9SH London  
United Kingdom

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**smilesonic**

SONIC TOOTHBRUSH

Thank you for purchasing the **smilesonic** sonic toothbrush. Before using the product, make sure to read the instruction manual carefully and save it in case you are going to need it later. The manufacturer does not hold responsibility for improper maintenance and incorrect use of the device.

## Specifications and contents of the set

Product type:	<b>smilesonic EX</b>
Working modes:	Clean/White/Polish/ Gum Care/Sensitive
Switch type:	A single button
Total working time:	60-180 minutes
Charging time:	9-16 hours
Operating voltage:	4V
Rated power:	3W
Charging current:	40-70 mA
Battery:	Lithium/14500/600 mAh
Battery indicator:	Flashing light
Waterproof rating:	IPX7

### The set contains:

Toothbrush handle, 4 brush heads, travel case, charger.

Protective cap

Brush head

Toothbrush handle

Power/mode change button

A panel displaying currently chosen mode

Battery indicator

Stand with a built-in inductive charger

## Brush head types

### DAILYCARE

Replaceable brush heads for complex care of teeth and entire oral cavity. Full protection and fresh breath guaranteed every day.

### SENSITOUCH

Replaceable brush heads for gentle brushing of sensitive teeth and gums. Full comfort and maximum effectivity of brushing guaranteed.

### SHINYWHITE

Replaceable brush heads for restoring the natural white color of the teeth. Effective at removing both the residue and discolorations.

### ORTHOCLEAN


Replaceable brush heads for advanced cleaning of interdental surfaces and places hardest to reach. Highly effective for cleaning dental braces and implants.

# Usage


## 01 Assembly

Before using the device, mount the brush head on top of the handle. There should be a small gap between the handle and the head. Avoid applying too much force while attaching the head to the handle.
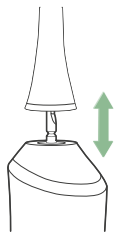
## 02 Rinsing

Before every brushing, rinse the bristle with running water and apply the right amount of toothpaste on it.
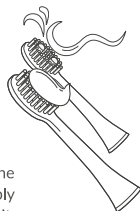
## 03 Choice

The toothbrush offers five cleaning modes, allowing you to adjust it to your needs.

**Tip:** It is recommended to choose the Sensitive mode during your first brushing.

### Clean

62000/min  
Removes tartar  
and bacteria

### White

96000/min  
Tooth whitening,  
removal of stains  
and discoloration

### Polish

70000/min  
Tooth polishing  
and brightening

### Gum Care

62000/min  
Massages the gums  
and improves their  
blood circulation


### Sensitive

62000/min  
Gentle brushing of  
sensitive teeth  
and surfaces

## 04 Brushing


The bristle should always be inclined at a 45 degree angle, in order to clean every spot within our mouth cavity.

**Tip:** It is recommended to turn the device on after putting it in our mouth first.
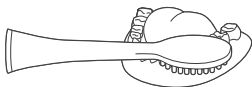
## 05 Rinsing

Wash your mouth out with water, then do the same with the toothbrush, removing bacteria and the remaining toothpaste.
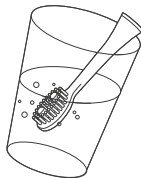
## Properties

### Remembering the last mode chosen

After turning the device on again, it will be set to the previously chosen cleaning mode.

### Two-minute timer

The toothbrush can turn itself off automatically. This helps us to brush our teeth for two minutes, just like the dentists recommend.

### Thirty-second interval

The toothbrush briefly stops working every 30 seconds, in order to let us know how much time has passed. This allows us to keep track of brushing every quarter of our oral cavity every 30 seconds.

### Indicator bristles

When the bristle starts losing its color, it means that the brush head is getting worn out. It is recommended to replace the brush head every three months.

### Up to 180 minutes of work

After a single charge, the toothbrush can be used for two months (work in Sensitive mode).

### One button

The device is operated by a single button. Pressing the button continuously after turning the toothbrush on allows us to change the cleaning mode.

## Correct brushing

1. Slowly and gently move the brush head along the teeth in a back-and-forth motion, allowing the bristle to get to the interdental surfaces.
2. Do not brush the teeth the same way you did while using the manual toothbrush.
3. In order to increase the effectivity of brushing, gently apply the toothbrush to the teeth, allowing it to brush them for you.
4. It is recommended to brush not only your teeth, but gums as well. The suggested order of toothbrushing is: the outer surface, the inner surface and finally the chewing surface.

## Cleaning the toothbrush

1. Clean the toothbrush regularly. Wash the handle and brush head with running water after each use.
2. Carefully dry all parts before putting them on the stand.
3. Before cleaning the charger, make sure to unplug it.
4. It is forbidden to clean the charger under the running water or soaking it.
5. Carefully dry the charger before plugging it back on.

## Charging

1. Before using the toothbrush for the first time, make sure that the device is fully charged.
2. Plug in the charger and then put the toothbrush in it.
3. During the charging procedure, the battery indicator should blink.
4. The indicator will continue blinking until the toothbrush is fully charged.
5. Once charged, the device will indicate the full battery by flashing light.
6. A full charging cycle lasts from 9 to 16 hours, allowing for 60 to 180 minutes of continuous work. (depending on the selected brushing mode)
7. We can unplug the charger or leave it be (the device is protected against overcharging).
8. Rapid blinking of the battery indicator means that the battery power is running low.

## Safety

In order to avoid damages or injuries, it is important to follow the basic safety rules related to using the electronic devices, including the following:

1. Before using the device for the very first time, make sure to read its instruction manual.
2. Before plugging the device, make sure that the voltage shown on the device is the same as the one in the outlet.
3. If the device is damaged in any way or makes unusual noises after being powered on (brush head, handle, charger, travel case), do not use it. It is forbidden to try to repair it on your own. The device does not contain any elements made for the user's repairs.
4. If the charger, travel case or brush heads are damaged, replace them with the original parts of the same type, in order to avoid any potential danger.
5. While unplugging the device, always grab and pull the plug, not the cable.

6. The charger and power cord should be used and stored away from the heat sources, water and damp.
7. The device should not be used by the people (including children) with physical, sensory or mental impairments and those lacking knowledge or experience necessary to operate such devices. The children should be instructed to not use the device as a toy.
8. The device should not be used in any way other than the one specified in the manual. The toothbrush is made only for cleaning teeth, gums and tongue.
9. The manufacturer does not hold responsibility for the damages caused by an improper use of the product.
10. Do not use bootleg accessories.
11. Do not use the toothbrush and the charger while in bath or under the shower.
12. Do not use this device in places where it can land or be thrown into the water.
13. Do not submerge the handle or the charger in water or any other liquid. After cleaning the device, make sure that every component is dry before plugging it back.

14. The device should not be used with wet hands.
15. The device is for the home use only. It is not made for the professional use.
16. If you had to undergo a jaw/gum surgery within the last 2 months, make sure to contact your dentist before using the device.
17. Consult your dentist if your gums start bleeding after using the toothbrush.
18. If you are using a pacemaker or a similar device, contact your dentist or the product's manufacturer before using the toothbrush.
19. Replace the brush head every 3 months or sooner, if you will notice the signs of wearing out.
20. Do not clean the brush heads in the dishwasher.
21. Do not use the brush heads other than the ones recommended by the manufacturer.

**CAUTION:** The device has a waterproof handle and is secured against the damp, but it should never be used underwater!

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Correct disposal of this product (waste electrical & electronic equipment)

After the product has reached the end of its service life, it is forbidden to throw it away along with the regular waste. Instead, it should be brought to your local collection/recycling point for the electric and electronic devices. It is indicated by the symbol found on the product, packaging or in the manual. The materials used for manufacturing the device can be recycled and used again. Thanks to recycling, reusing or other forms of utilizing the products past their service life, you are all contributing to protecting the environment. The information about finding the correct recycling point can be obtained at your local municipal administration or retailer we purchased the product at.

## smilesonic Ltd

C220d Trident Business Centre  
89 Bickersteth Road  
SW17 9SH London  
United Kingdom  
info@smilesonic.co.uk

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